


U-Tube Type Min-Max Thermometer

This thermometer can measure minimum, maximum and current temperatures.


Features

- The min. and max. temperatures can be read with the indexes.
- The min. and max. temperature readings are easily reset by pushing reset button.

Reading temperatures

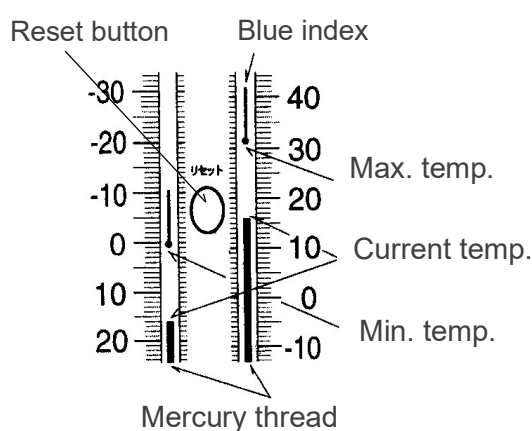

- Max. temp. : lower end of index on the right
- Min. temp. : lower end of index on the left
- Current temp. : both upper ends of the mercury threads

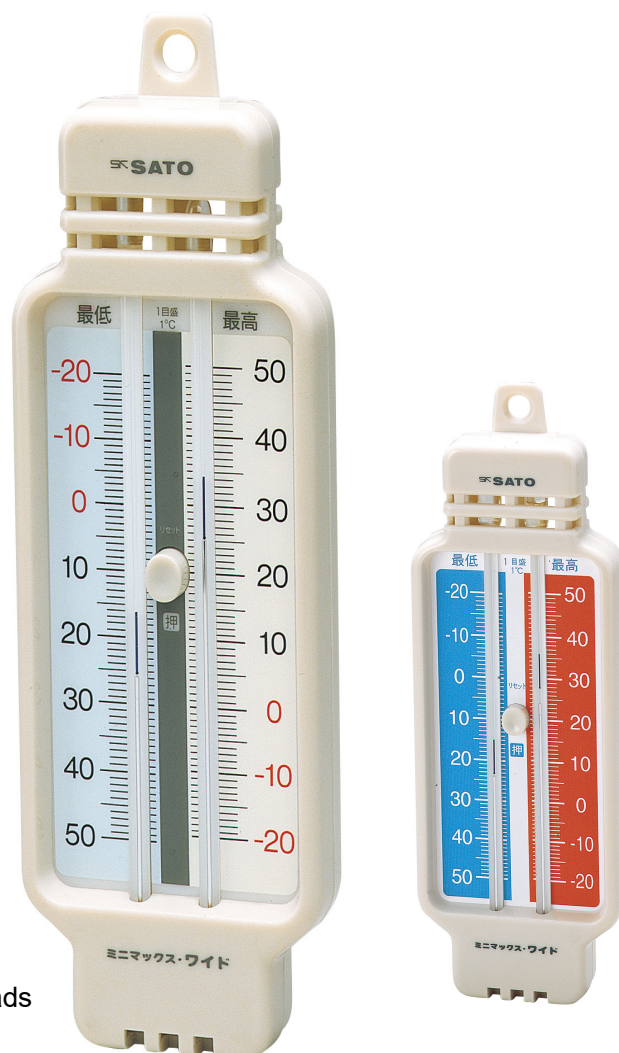

Specifications

Cat. No.	1506-00	1506-50	1508-00	1508-50
Scale color	Ivory for max. temperature Light-blue for min. temperature		Red for max. temperature Blue for min. temperature	
Measuring range	-20 to 50°C	-40 to 50°C	-20 to 50°C	-40 to 50°C
Min. graduation	1°C			
Sensing liquid	Mercury, Creosote			
Materials	Thermometer : Glass Body : ABS resin			
Dimensions	(W)76 x (H)220 x (D)25mm			
Weights	approx.135g			

SATO KEIRYOKI MFG. CO., LTD.

3-4, Kanda-kajicho, Chiyoda-ku, Tokyo 101-0045 Japan

URL: <http://www.sksato.co.jp/english/>